


RiMed Inc, your reliable connection for AFT beds, introduces the

# HydroAire

GENERATION II

## Air-Fluidized Therapy (AFT) Bed System

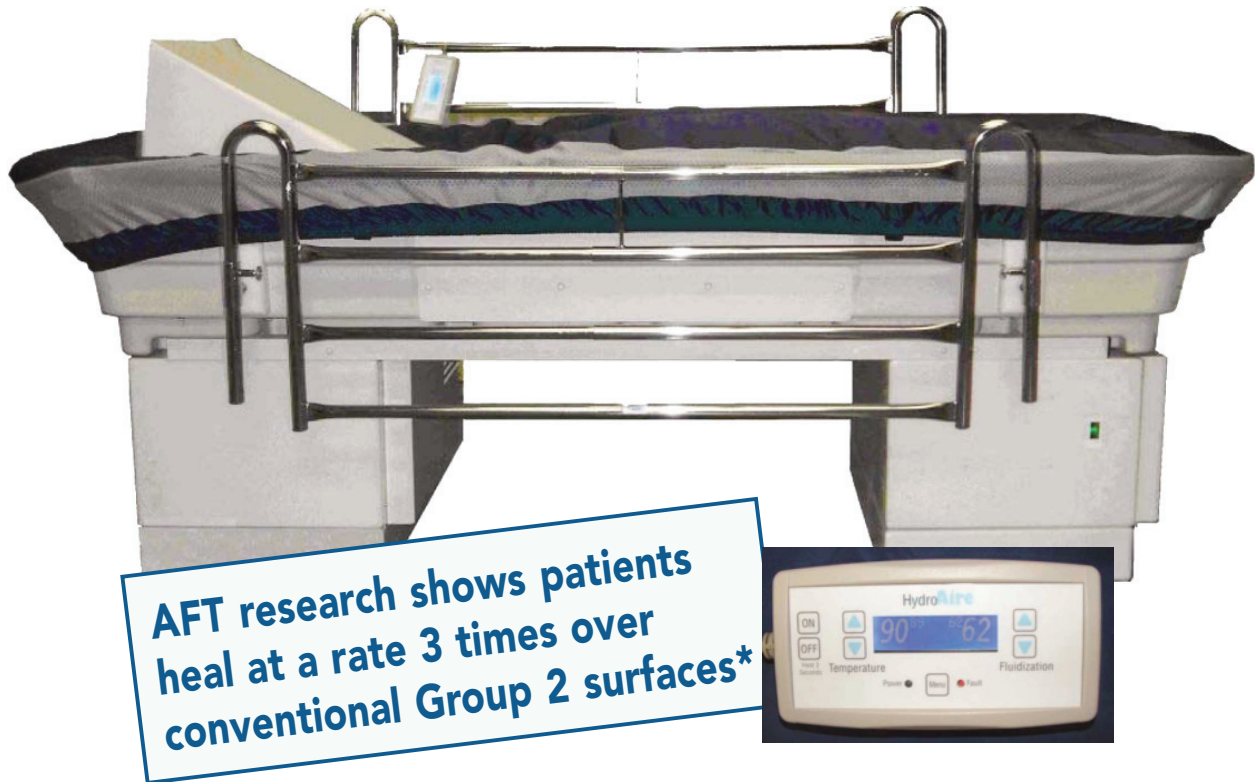

AFT research shows patients heal at a rate 3 times over conventional Group 2 surfaces\*



**RiMed Inc.**

**Your  
reliable connection  
for AFT beds!**

**Call today for more  
information: 1-866-942-2650  
www.rimedinc.com**

- AFT – the Gold Standard in surfaces for patients with non-healing Stage 3 or 4 pressure ulcers on the trunk or pelvis – is a proven, reliable solution for home care.
- Utilization of AFT results in better outcomes and allows agencies to reduce visits and increase healing rates.
- The HydroAire represents the best in AFT: Better fluidization, quieter, creates less heat – the **entire surface fluidizes**. It is suitable for mobile homes. *Learn more about the HydroAire AFT at [www.auroramfg.com](http://www.auroramfg.com)*
- Reduces electricity costs by 40% over previous unit. RiMed Inc. is a consistent provider of AFT to Medicare patients.

\*Comparison of air-fluidized therapy with other support surfaces used to treat pressure ulcers in nursing home residents. Dr. Ochs, et al. *Ostomy Wound Management*, Feb. 2005.

# HydroAire

GENERATION II

## AFT: The best surface for non-healing Stage 3 or 4 pressure ulcers

- Pay For Performance! AFT is a solution for chronic pressure wounds.
- AFT is best solution for patients who have not improved on Group 2 surfaces, and a potential solution for patients 'capping out' on Group 2 surfaces.
- Patients who have had myocutaneous skin flap procedures are great candidates!
- Reduces risk of acquired ulcers.
- Reduces incidence of re-hospitalization.\*
- RiMed consultants visit with patients before and after delivery to ensure compliance
- RiMed provides same-day response for servicing beds.

**Call today for more information: 1-866-942-2650**  
**www.rimedinc.com**



### How It Works:

\*Patients body envelops into bead media, beads conform to the patients body, evenly distributing their weight.

\*\*With pressures under 32mmg it eliminates the need for repositioning.

\*\*\*Provides the best environment to enhance blood flow and wound healing

### Product Features for Model 3000

- Adjustable height enhances patient mobility
- Modular design for easy in-home installation
- Breathable GoreTex® sheets
- Dual blower design provides redundancy in event of blower failure
- Adjustable Fluidization
- 84" long x 39" wide; Height adjustable 26" -30"
- Variable heating temp 80° to 100°
- Remote Control to change parameters
- Bed Weight - empty 282 lbs.; w/ beads 1,122 lbs.
- Siderails are standard for safety
- Trapeze bar and stand available
- Works in conjunction with standard patient lifter
- Slant wedge available at varying angles
- Accessible blower assembly allows easier servicing

### Specifications

- Maximum recommended patient weight limit ..... 250 lbs.
- Bed Height ..... 30 or 26 inches
- Bed Width  
Overall ..... 36 inches  
Inside Tank ..... 33-3/4 inches
- Bed Length  
Overall ..... 84 inches  
Inside Tank ..... 81-3/4 inches
- Bed Weight  
Without Micro Spheres ..... 282 lbs.  
Micro Spheres ..... 840 lbs.  
With Micro Spheres ..... 1,122 lbs.  
Max Patient Size ..... 250 lbs.
- Electrical  
Voltage ..... 115 VAC  
Frequency ..... 60 Hertz  
Noise Level ..... <50 DBa  
Power Cord Length ..... 10 feet

NOTE: Federal law restricts this device to sale or rental by or on the order of a physician